



Knowledge Beyond Measure.

Aerosol Neutralizers

Models 3012, 3012A, 3054, 3054A, 3077, 3077A, 348002



Minimize particle losses and coagulation caused by electrostatic charges, or charge particles properly for scanning mobility particle sizer analysis and air-filter measurements

Aerosol particles dispersed by nebulization, combustion, or powder dispersion are usually electrostatically charged. Normally, a high level of electrical charge is undesirable. It increases the particle loss to the walls of transport and sampling systems, or it can affect filter-efficiency measurements. To ensure that instruments operating on the electrostatic principle work properly, aerosol particles entering these instruments must be neutralized.

Real-time viable particle detection enables

- Models 3054 and 3054A are designed to fit onto a TSI® Vibrating Orifice Aerosol Generator Model 3450
- Models 3077 and 3077A are designed to fit into TSI® series Electrostatic Classifiers Model 3080 and all Scanning Mobility Particle Sizer spectrometers
- Model 348002 is designed to fit into a TSI® Electro Spray Aerosol Generator Model 3480
- TSI Aerosol Neutralizers may be used for general-purpose neutralization. If your aerosol carries higher charge levels or if your aerosol flow rate is in the upper half of the specified range, TSI recommends the 3012A, 3054A, or 3077A. See flow rate ranges and other specifications on the reverse side of this document


Features and Benefits

- Neutralize charge on aerosol particles
- Increase measurement efficiency
 - Displays up to three environmental parameters
 - Stainless steel enclosure

Operation

TSI Aerosol Neutralizers use a radioactive source (Kr-85 or Po-210) to perform this function. The radioactive source ionizes the surrounding atmosphere into positive and negative ions. Particles carrying a high charge can discharge by capturing ions of opposite polarity. After a short time, the particles reach charge equilibrium such that the aerosol carries a bipolar distribution as described in the literature by Fuchs¹.

In the longer, cylindrical neutralizers (shown on first page), inert Kr-85 gas is sealed completely inside a stainless-steel tube with wall thickness of 0.05 millimeters. The tube is shielded by a metal outer housing. As described by Liu and Pui², aerosol passes through the space between the stainless-steel tube and the outer housing. Regarding the Model 348002 shown at right, aerosol does not pass through the neutralizer. Instead, the Model 348002 is placed adjacent to an aerosol stream.

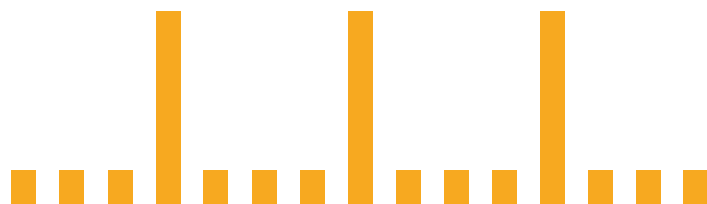

Specifications

Aerosol Neutralizers

Models 3012, 3012A, 3054, 3054A, 3077, 3077A, 348002

Model	3012	3012A	3054	3054A	3077	3077A	348002
Source	Kr-85						
Emission	Beta decay						
Radioactivity [MBq (mCi)]	74 (2)	370 (10)	370 (10)	740 (20)	74 (2)	370 (10)	185 (5)
Half-life	10.7 years						
Maximum Flow Rate	50 liters/min		150 liters/min		5 liters/min		
Maximum Temperature	50°C (122°F)						
Maximum Pressure	35 kPa (5 psig)						
Outer Housing Material	Anodized aluminum		Stainless steel			Aluminum	
Dimensions							
Inlet Diameter	1.27 cm (0.5 in.)		2.46 cm (0.97 in.)		0.64 cm (0.25 in.)		NA
Outlet Diameter	1.27 cm (0.5 in.)		3.50 cm (1.38 in.)		0.64 cm (0.25 in.)		NA
Housing Diameter	7.72 cm (3.04 in.)		8.94 cm (3.52 in.)		3.89 cm (1.53 in.)		3.20 cm (1.26 in.)
Overall Length	52.86 cm (20.81 in.)		64.14 cm (25.25 in.)		21.44 cm (8.44 in.)		0.97 cm (0.38 in.)
Weight	1.0 kg (2.2 lb)		3.5 kg (7.7 lb)		0.4 kg (0.9 lb)		14 g (0.5 oz)

Relocation and Disposal

Aerosol Neutralizers contain a radioactive source that requires special handling. Government regulations require that you notify TSI® Incorporated if you wish to dispose of this product or transfer it to a different facility. TSI® is required by law to assist you with this process. There is a fee for this service. Please contact TSI® customer service for additional information.

Bibliography

¹Fuchs NA, "On the Stationary Charge Distribution on Aerosol Particles in Bipolar Ionic Atmosphere," Geophys. Pura Appl. 56:185 (1963).

²Liu BYH, and DYH Pui, "Electrical Neutralization of Aerosols," J. Aerosol Sci. 5:465-472 (1974).

TSI has been issued license number 1154-200-62 by the Minnesota Department of Health to sell and distribute these Aerosol Neutralizers. Users in the United States need not apply for additional U.S. Government licenses to handle these products. However, some state and local governments may require special licenses and some organizations may have special handling procedures. Check all local requirements.

To Order Specify

3012

3012A

3054

3054A

3077

3077A

340082

Description

For general-purpose aerosol applications with high flow rates; has screened inlet and outlet

For general-purpose aerosol applications with high flow rates, has screened inlet and outlet and five times the activity of Model 3012

For use with Model 3450 Vibrating Orifice Aerosol Generators

For use with Model 3450 Vibrating Orifice Aerosol Generators; has twice the activity of Model 3054

For general-purpose applications involving moderately charged aerosols or low flow rates; lower activity reduces regulatory concerns compared to Model 3077A

For general-purpose applications involving highly charged aerosols and/or high flow rates; used with Model 3034 Scanning Mobility Particle Sizer and other TSI instruments; has five times the activity of Model 3077

For use with Model 3480 Electro Spray Aerosol Generators

Specifications are subject to change without notice.

TSI and the TSI logo are registered trademarks of TSI Incorporated in the United States and may be protected under other country's trademark registrations.



Knowledge Beyond Measure.

TSI Incorporated - Visit our website www.tsi.com for more information.

USA Tel: +1 800 874 2811 India Tel: +91 80 67877200
 UK Tel: +44 149 4 459200 China Tel: +86 10 8219 7688
 France Tel: +33 1 41 19 21 99 Singapore Tel: +65 6595 6388
 Germany Tel: +49 241 523030