




Knowledge Beyond Measure.

Qualitative Fit Testing with FitPro™ Ultra Fit Test Software

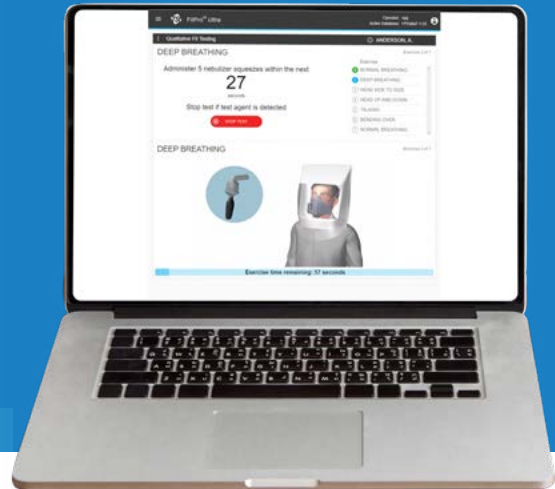

Your Main Screen to Qualitative Fit Testing

Enhance your fit testing process with step-by-step guidance using Qualitative Fit Testing with Software FitPro™ Ultra Fit Test Software from TSI® to help ensure you're doing qualitative fit testing consistently and correctly.

Consistent Testing Protocol

- Visually guided OSHA 7-step exercise protocol for a proper and complete qualitative fit test
- Easy-to-follow on-screen instructions, test agent delivery timers with exercise fit test animations

Intuitive Administrator Guidance

- Automatically guides the fit test administrator step-by-step through the fit test process
- Indicators for proper test solution dispensing to deliver the right number of squeezes at the right times

Automated Record Keeping and Reporting

- Storage of fit test results in a database for easy retrieval for printing fit test records and report generation

Additional Features and Benefits

- Automates fit test process using built-in fit test exercise
- Local file based default database (optional support for server-based Microsoft® SQL Server®)
- Comprehensive database management functionality
- Maintains separate databases for different companies/departments/locations
- Wide-range of reporting and printing capabilities, including individual fit test reports, due or past due fit test reports, custom reports, and printing of individualized fit test cards (using card-form paper)

Single Software Solution for Both Fit Test Methods

Built on a proven platform, experience the simplicity the FitPro™ Ultra Fit Test Software for Qualitative Fit Testing can bring to your respiratory fit test program. For those also completing quantitative fit testing with PortaCount® Respirator Fit Tester you now can add the qualitative fit testing option using the same FitPro™ Ultra Fit Testing Software. It provides a seamless fit test administrator experience and consolidated recording keeping in a single software solution.



Specifications

Qualitative Fit Testing with FitPro™ Ultra Fit Test Software

Computer Requirements

- Microsoft® Windows® 10 (32-bit & 64-bit) or Windows 11 operating systems
- Microsoft .Net version 4 (installation software will install Microsoft.NET automatically, if necessary)
- Note—An Internet connection is required to install Microsoft.NET
- >1 GHz processor
- >1280x800 screen resolution monitor
- 500 MB free space on hard drive
- Internet connection (to download the software)

Specifications are subject to change without notice.

Licensed under US Patent No. 5,117,190 and other patents from Commissariat a L'Energie Atomique, France. U.S. Patent No. 6,125,845. Other patents pending. Access, Microsoft, SQL Server and Windows are either registered trademarks or trademarks of Microsoft Corporation in the United States and/or other countries.

PortaCount, TSI and the TSI logo are registered trademarks in the United States and may be protected under other country's trademark registrations.



Knowledge Beyond Measure.

TSI Incorporated - Visit our website www.tsi.com for more information.

USA	Tel: +1 800 874 2811	India	Tel: +91 80 67877200
UK	Tel: +44 149 4 459200	China	Tel: +86 10 8219 7688
France	Tel: +33 1 41 19 21 99	Singapore	Tel: +65 6595 6388
Germany	Tel: +49 241 523030		