

FITPRO+™ FIT TEST SOFTWARE FOR PORTACOUNT® PRO AND PRO+ RESPIRATOR FIT TESTERS

FitPro+™ Fit Test Software makes fit testing simple and foolproof. Included with every PortaCount Pro and PortaCount Pro+ Respirator Fit Tester, the program automatically leads you through the fit test protocol step by step – you can't make a mistake. Your fit test results are stored in a database for easy report generation, retrieval and regulatory compliance.



Computer Requirements


- + Microsoft® Windows® 7 (32 and 64 bit), Windows 8 (32 and 64 bit), or Windows 10 (32 and 64 bit) operating system
- + Microsoft .Net version 4 (installation software will install Microsoft .NET automatically, if necessary)
Note—An Internet connection is required to install Microsoft .NET
- + >1 GHz processor
- + >1280x800 screen resolution monitor
- + 500 MB free space on hard drive
- + 1 free USB port
- + CD-ROM reader or internet connection (to install the software)

Features and Benefits

- + Operate up to four (4) PortaCounts simultaneously using one computer (PortaCount Models 8030/8038 only)
- + Select between new Horizontal or Vertical Fit Test Dashboard
- + Automates fit test process using built-in fit test exercise protocols
- + Microsoft Access® default database (optional support for server-based Microsoft SQL Server®)
- + User Interface in five languages (English, German, Polish, French and Spanish)
- + Comprehensive database management functionality (merging, moving, copying and deleting records; exporting databases; and backup and restoring of data)
- + Maintains separate databases for different departments/locations
- + Easily transfers records between PortaCount Fit Testers using a flash drive
- + Wide-range of reporting and printing capabilities, including individual fit test reports, due or past due fit test reports, custom reports, and printing of individualized fit test cards (using card-form paper or optional card printer)
- + Real-time fit factor display and fit testing emulators for respirator training and familiarization



UNDERSTANDING, ACCELERATED

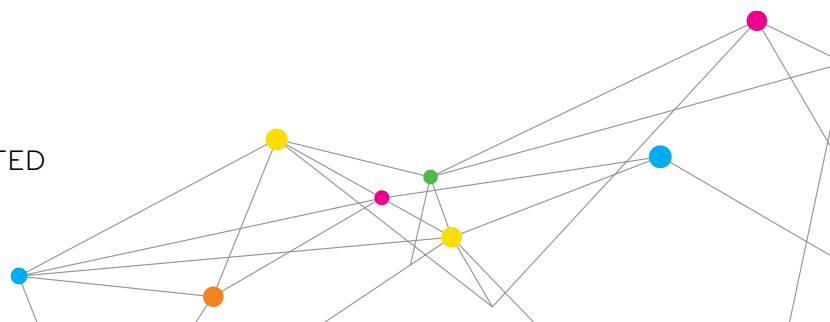

SPECIFICATIONS

FITPRO+ FIT TEST SOFTWARE FOR PORTACOUNT PRO AND PRO+ RESPIRATOR FIT TESTERS

Fit Test Report and Card Printing Options

FitPro+ is designed to improve your respirator program and help protect workers through easy generation of individualized fit test reports inclusive of personal wallet cards. The wallet cards provide documentation of the completed fit test and easy field verification of proper respirator size and model for authorized respirator users.

The reports can be printed by plugging a USB printer directly into the PortaCount Pro/Pro+ Respirator Fit Tester. Exclusive card-form paper is available direct from TSI that provides for the convenient, durable tear-out wallet card.

FitPro+ Fit Test Software also maintains the ability to print fit test results to an external ID card printer.

FIT TEST REPORT			
3132014	ID NUMBER	1234	
	LAST NAME	BRUCE	
	FIRST NAME	MORTON	
	COMPANY		
	LOCATION		
	TEST DATE	1/22/2014 10:52 AM	PORTACOUNT SN
	DUPLICATE	1/22/2015	003005502
			NON COMPANION
			N
	RESPIRATOR	MSA ULTRATWIN FULL FACE (500)	PROTOCOL
	MANUFACTURER	MSA	OSHA 29CFR1910.134
	MODEL	ULTRATWIN	PASS LEVEL
	MASK STYLE	FULL FACE	500
	MASK SIZE	F	APPROVAL
			EFFICIENCY=95%
			Pass
			False
	EXERCISE	DURATION(msec)	FIT FACTOR
	NORMAL BREATHING	64	110991
	DEEP BREATHING	64	34756
	HEAD SIDE TO SIDE	64	219443
	HEAD UP AND DOWN	64	63325
	TALKING	64	70031
	GRINACE	15	6165
	BENDING OVER	64	165212
	NORMAL BREATHING	64	64501
	OVERALL FF		77582
			Y
	FITTEST OPERATOR	JDM	DATE
	NAME	MORTON BRUCE	DATE
	NOTE:		

Respirator Fit Test Card


Name: MORTON BRUCE Test Date: 1/22/2014
 ID: 1234 Next Test Date: 1/22/2015

Respirator	Results
Mfg: MSA	Overall FF: 77582
Model: ULTRATWIN	FF Pass Level: 500
Style: FULL FACE	Pass: Y
Size: F	Operator: JDM
	Protocol: OSHA 29CFR1910.134
	Fit Test Method: QNFT using TSI PortaCount
	*** Your company contact information here ***

Licensed under US Patent No. 5,117,190 and other patents from Commissariat a L'Energie Atomique, France. U.S. Patent No. 6,125,845. Other patents pending.

PortaCount, TSI and the TSI logo are registered trademarks, and FitPro is a trademark of TSI Incorporated.

Access, Microsoft, SQL Server and Windows are either registered trademarks or trademarks of Microsoft Corporation in the United States and/or other countries.

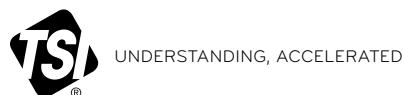

TSI Incorporated - Visit our website www.tsi.com for more information.

USA	Tel: +1 800 874 2811	India	Tel: +91 80 67877200
UK	Tel: +44 149 4 459200	China	Tel: +86 10 8219 7688
France	Tel: +33 1 41 19 21 99	Singapore	Tel: +65 6595 6388
Germany	Tel: +49 241 523030		